

Northwest Life Science Specialties, LLC

Premier Products for Superior Life Science Research

NWLSSTTM *Urinary 8-OHdG* *ELISA*

Product NWK-80HDG01
For Research Use Only



Simple assay kit for measurement of 8-hydroxydeoxyguanosine (8-OHDG) in urine samples.

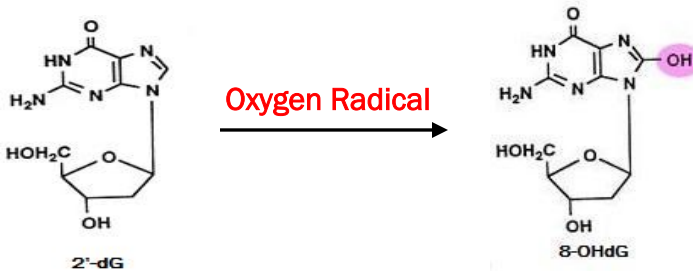
Samples such as plasma, tissue homogenates and cell lysates typically contain lower amounts of 8-OHDG than urine; we recommend using the more sensitive NWK-80HDG02 assay for these sample types.

Table of Contents

Section	Page
Introduction	3
Intended Use	3
Test Principle	3
General Specifications	4
Kit Contents	4
Required Materials Not Provided	4
Required Instrumentation	4
Warnings, Limitations, Precautions	5
Storage Instructions	5
Assay Preparation	6
Reagent Preparation	6
Sample Handling/Preparation	6
Assay Protocol	7
Data Analysis	8
References	8
Statement of Limited Warranty	9
End-User Notes	10

Introduction:

Oxidative stress is known to play an important role in the development of various diseases and aging process. Under conditions of oxidative stress, 8-hydroxy-2'-deoxyguanosine (8-OHdG) is formed when DNA is oxidatively damaged by reactive oxygen species (ROS). 8-OHdG is one of the most sensitive biomarkers for oxidative stress and can be detected in serum, urine and/or DNA isolated from cells and tissues in humans and animals. The figure below depicts the oxidation 2'-deoxyguanosine (2'dG) to form the 8-OHdG biomarker.

**Intended Use:**

The NWLSS™ 8-OHdG ELISA is intended for quantitative detection of the 8-hydroxyl 2' deoxyguanosine oxidative DNA adduct in biological samples including urine, plasma, serum, tissue and other biological samples.

Test Principle:

The NWLSS™ 8-OHdG ELISA assay uses a competitive format wherein a murine monoclonal antibody to 8-OHdG (*Primary Anti-8OHdG*) and sample or standard are added to a microtiter plate which has been precoated with 8-OHdG. Sample or calibrator 8-OHdG competes with plate-bound 8-OHdG for binding with the antibody. Accordingly, higher concentrations of sample or calibrator leads to reduced binding of the antibody to the 8OHdG coated plate. A subsequent wash step removes any free 8-OHdG/antibody adduct leaving stationary plate bound 8-OHdG complexed to antibody for later detection. Anti-murine antibody conjugated to horse radish peroxidase (*HRP-Conjugate*) is then added to the plate. HRP-conjugate binds to remaining murine anti-8-OHdG and unbound HRP-conjugate is removed in another wash step. Addition of 3,3',5,5'-tetramethylbenzidine (*TMB Substrate*) results in blue color development proportional to the amount of anti 8OHdG antibody bound to the plate and inversely proportional to the concentration 8-OHdG in original samples or calibrators applied to the plate. The reaction is terminated by addition of phosphoric acid (*Stop Solution*) producing yellow color with measurable absorbance at 450 nm.

General Specifications:

Format::	1 X 96 wells	
Number of tests:	Triplicate =	24
	Duplicate =	40
Specificity:	8-Hydroxy 2' deoxyguanosine	
Sensitivity:	LLD = 0.5 ng/mL in sample applied to plate	

Kit Contents:

Microtiter Plate:	Precoated with 8OHdG	12 X 8 wells
8OHdG Standard:	6 levels Purified 8OHdG	6 Vials
Primary Antibody:	Mu-Anti-8OHdG Monoclonal	1 vial (Lyoph.)
Primary Antibody Buffer:	Phosphate buffered saline	1 vial (6mL)
Secondary Antibody:	Anti-murine-HRP conjugate	1 vial (Lyoph.)
Secondary Antibody Buffer:	Phosphate buffered saline	1 vial (12mL)
TMB Substrate:	3,3',5,5'-tetramethylbenzidine	1 vial (0.25mL)
Diluting Buffer:	H ₂ O ₂ /Citrate/PBS	1 vial (12mL)
Wash Buffer (5x):	Concentrated PBS	2 vials (26mL)
Stop Solution:	1M Phosphoric acid	1 vial (12mL)
Plate Seal		2 Sheets

Required Materials Not Provided:

Adjustable pipettes capable of transferring 50 µL to 200 µL volumes.
 A multi-channel or repeater pipette (recommended).
 Reagent Trays (if using multi-channel).
 37 °C Incubator
 Distilled water
 PBS (pH 7.4)

Required Instrumentation:

Microplate reader with 450nm capability.

Warnings, Limitations, Precautions:

Incubation Temperature

Measured values can be affected by variations in incubation temperatures. Care should be exercised to assure uniform temperature during incubations, particularly during the Primary Antibody reaction period. It is recommended that incubation be performed in a humid environment (i.e. in an incubator with water present).

Sample pH

The optimal solution pH of sample mixed with primary antibody should be between 6.0 to 8.0. Samples with unusual pH can be diluted in PBS if pH is deemed to be an issue.

Possible edge effects

To minimize edge effects, ensure that plate is sealed properly and that incubation times are as uniform as possible. To maintain the most uniform temperature within the wells, it is recommended that any unused wells on a single strip be filled with an equal volume of water prior to incubation.

TMB Substrate:

Do not dilute TMB Substrate until just before needed. Keep it in the dark.

Blank wells:

Do not add Primary Antibody to Blank wells.

Storage Instructions:

Upon receipt, store the reagents at 2-8°C. This kit is stable if unopened for 9 months from the date of manufacture. Do not use components beyond the expiration date printed on the kit box label.

Unused wells and opened reagents must be stored at 2-8°C and should be used within two weeks after first opening.

Assay Preparation*Plate Reader Setup*

Wavelength:	450 nm
Mode:	Endpoint
Shaker:	On
Curve Fit:	4-Parameter Logistic

Reagent Preparation:

Bring all reagents, plate wells to be used samples and calibrators to room temperature (20-25 °C) before use.

Wash Buffer:

Prepare the necessary volume of wash buffer by mixing 1 part **5X Wash Buffer** with 4 parts distilled water. Label as **Working Wash Buffer**.

Primary Antibody: (Red Cap and Label)

Reconstitute the **Primary Antibody** with **Primary Antibody Buffer**.

Secondary Antibody: (Blue Cap and Label)

Reconstitute the **Secondary Antibody** with **Secondary Antibody Buffer**.

TMB Substrate:

Just prior to use, prepare the necessary volume of TMB Substrate by mixing 1 part **TMB Substrate** with 100 parts **Dilution Buffer**. Label as **Working TMB Substrate**.

Sample Handling/Preparation:

The multi-disciplinary interest in measuring 8-OHdG has resulted in a myriad of sample types and experimental conditions and is beyond the scope of this product insert to describe sample processing in detail for each case. However, general guidelines are provided below for representative sample types.

Urine:


If samples are cloudy (opaque) they should be centrifuged at ~2000-5,000g for 10-15 minutes. If urine is clear, no pretreatment is necessary and samples may be applied directly to the assay. "Abnormal" urine may require a 3X dilution with PBS (pH 7.4) prior to assay to ensure an assay pH between 6.0 and 8.0.

Assay Protocol:*Standard Procedure for Microplate Assay*

1. Bring all reagents to room temperature.
2. Remove appropriate number of wells for assay from foil pouch.
3. Construct an assay template to ensure proper sample addition.
4. Add 50 μ L of **Sample** or **Standard** to each well to be assayed.
5. Add 50 μ L **PBS** to blank wells.
6. Add 50 μ L **Reconstituted Primary Antibody** to each well to be assayed except blank wells.
7. Shake lightly side to side to ensure proper mixing.
8. Cover plate with adhesive strip, then incubate at 37 °C for 1 hour.
9. Empty contents of wells into sink and blot on paper towel to remove as much fluid as possible.
10. Wash plate 3 times as follows:
Add 250 μ L **Working Wash Buffer** per well
Shake plate slightly during soak period for best wash results.
Empty wash solution into sink by inversion then blot plate against clean paper towel to remove any remaining washing buffer.
11. Add 100 μ L of **Reconstituted Secondary Antibody** per well.
12. Shake lightly side to side to ensure to mixing.
13. Cover the plate with an adhesive strip then incubate at 37 °C for 1 hour.
14. Empty contents of wells into sink and blot on paper towel.
15. Wash the plate 3 times as in Step 10.
16. Add 100 μ L of **Working TMB Substrate** per well.
17. Shake lightly side to side to ensure proper mixing.
18. Incubate for 15 minutes, at room temperature in the dark.
19. Add 100 μ L of **Stop Solution** per well.
20. Measure the absorbance at 450nm.

Data Analysis:

Create a standard curve by plotting Absorbance vs. Concentration (log scale) for each standard level assayed. If available, set the plate reader to utilize 4-Parameter curve fit. An example standard curve is shown below.

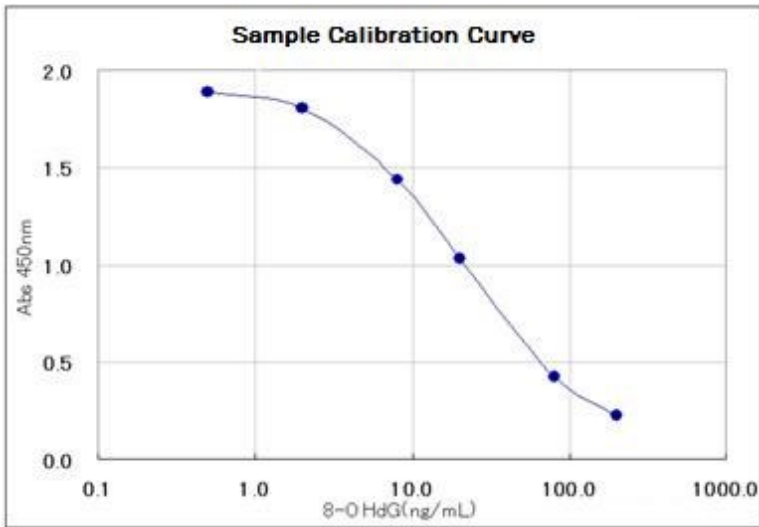**References:**

1. S.Okamoto, and H.Ochi (1992) *Chemical Abst.*; 129859a.
2. H.Kasai, P.F.Crain, Y.Kuchino, S.Nishimura, A.Ootsuyama, and H.Tanooka (1986) *Carcinogenesis* 7; 18491851
3. S.Toyokuni, T.Tanaka, Y.Hattori, Y.Nishiyama, A.Yoshida, K.Uchida, H.Hiai, H.Ochi, and T.Osawa (1997) *Lab.Invest.* 76, 365374
4. M.D.Evans, M.S.Cooke, I.D.Podmore, Q.Zheng, K.E.Herbert, and J.Lunec (1999) *Biochemical and Biophysical Research Communications* 259 pp374378; Discrepancies in the measurement of UVCinduced 8-oxo 2' deoxyguanosine: Implications for the analysis of oxidative DNA damage.
5. Tomoko Shimoike, Toyoshi Inoguchi, Fumio Umeda, Hajime Nawata, Katsumi Kawano and Hiroto Ochi (2000) *Metabolism* 49(8) pp10301035;
The meaning of serum levels of advanced glycosylation end products in diabetic nephropathy.

References (Continued):

6. WenYing Fan, Kazunori Ogusu, Katsuyasu Kouda, Harunobu Nakamura, Tomoaki Satoh, Hirotoomo Ochi and Hiroichi Takeuchi (2000) *J Physiol Anthropol* 19(6) pp287289; Reduced oxidative DNA damage by vegetable juice intake: A controlled trial.

7. H.Ochi, M.Hashimoto, J.Kurashige (2001) *Culture and Science*; Proceedings of 2001 International Conference on OCha (tea).

8. M.S.Cooke, M.D.Evans and J.Lunec (2002) *Free Radical Research* 36(9) pp929932; DNA repair: insights from urinary lesion analysis.


9. Tadashi Matsubasa, Takako Uchino, Shinnyo Karashima, Yuichi Kondo, Kenichi Maruyama, Masako Tanimura and Fumio Endo (2002) *Free Radical Research* 36(2) pp189193; Oxidative Stress in very low Birth Weight Infants As Measured by Urinary 8OHdG.

Statement of Limited Warranty:

Northwest Life Science Specialties, LLC (NWLSS) makes no guarantee of any kind, expressed or implied, that extends beyond the description of the material in this kit, except that they will meet our specifications at the time of delivery. Customer's remedy and NWLSS' sole liability is limited to, at NWLSS' option, refund of the purchase price, or the replacement of material not meeting our specification. By acceptance of our product, customer assumes all liability and will indemnify and hold NWLSS harmless for the consequence of this product's use or misuse by the customer, its employees, or others. Refund or replacement is conditioned of customer notifying NWLSS within twenty-one (21) days of the receipt of product. Failure to give notice within 21 days shall constitute a waiver by the customer of all claims hereunder with respect to said product.

Notes:

Notes:

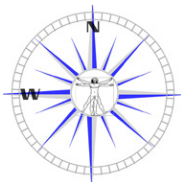

Northwest

Life Science Specialties, LLC

5131 NE 94th Avenue, Suite 201

Vancouver, WA 98662

Phone 360-449-3091 or Toll Free: 888-449-3091

Fax 360-449-3092

E-mail: sales@nwlifescience.com