


**Northwest**  
Life Science Specialties, LLC

## Anti Nitrotyrosine Polyclonal

<b>Catalog Number:</b>	NWA-NTRP0067
<b>Description:</b>	Rabbit anti Nitrotyrosine (100 µL per Vial).
<b>Specificity:</b>	Specific for Nitrotyrosine
<b>Format:</b>	100 uL IgG <sub>1K</sub> in HEPES with 0.15M NaCl, 0.01% BSA, 0.03% sodium azide, and 50% glycerol
<b>Antigen:</b>	Synthetic Peptide
<b>Reconstitution:</b>	NA
<b>Suggested Use:</b>	Western blotting (1:2,000) ELISA

User should determine optimum titer for  
Each application.

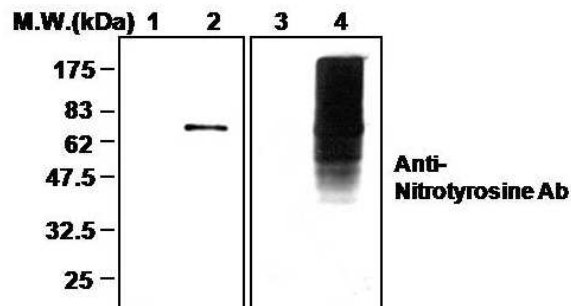

**Immunoblot Analysis of cell lysates**  
Lane 1 : Human albumin protein  
Lane 2 : Human albumin protein treated with 10mM  
peroxynitrite  
Lane 3 : A431 cell lysate  
Lane 4 : A431 cell lysate treated with 10mM peroxynitrite

<b>Cross Reactivity:</b>	Human, Rat, Mouse
<b>Storage:</b>	Store at -20 °C upon arrival.
<b>Stability:</b>	Store for 1 year at -20 degree from date of shipment
<b>References:</b>	1) Hanafy KA, et al, Med Sci Monit. 2001; vol.7(4): pp.801-19. 2) Nakazawa H, et al, Free Radic Res. 2000; vol.33(6): pp.771-84. 3) Ceriello A. Int J Clin Pract Suppl. 2002; vol.129: pp.51-8.

**For Research Use Only**

Northwest Life Science Specialties, LLC  
16420 S McGillivray, Suite 103, PMB 106, Vancouver, WA 98683  
Phone: (360) 449-3091, (888) 449-3091, Fax: (360) 449-3092, E-mail: [info@nwlifescience.com](mailto:info@nwlifescience.com)